Version 04/04

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© Copyright 2004 by Conrad Electronic GmbH. Printed in Germany.

100%  
recycling  
paper.  
Bleached  
without  
chlorine.

**GB Imprint**

These operating instructions are published by Conrad Electronic GmbH, Klaus-Conrad-Str. 1, 92240 Hirschau/Germany. No reproduction (including translation) is permitted in whole or part e.g. photocopy, microfilming or storage in electronic data processing equipment, without the express written consent of the publisher. The operating instructions reflect the current technical specifications at time of print. We reserve the right to change the technical or physical specifications.

© Copyright 2004 by Conrad Electronic GmbH. Printed in Germany.

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© Copyright 2004 by Conrad Electronic Ned BV. Printed in Germany.

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**(GB)** These operating instructions belong with this product. They contain important information for putting it into service and operating it. This should be noted also when this product is passed on to a third party.

Therefore look after these operating instructions for future reference!  
A list of contents with the corresponding page numbers can be found in the index on page 16.

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## Introduction

Dear customer,  
Thank you for buying the Body Fat Scales Bettina.  
With this device, you have purchased a product the design of which complies with the state-of-the-art.

**This product meets the requirements of the applicable European and national guidelines. Its conformity has been verified; the corresponding declarations and documents deposited with the manufacturer.**

To maintain this condition and to ensure risk-free operation, you, as the user, must observe these operating instructions”

***For a fast response of your technical enquiries please use the email address listed below.***

**Germany:**      **Fon. + 49 9604/40 88 80**  
**Fax + 49 9604/40 88 48**  
**e-mail: [tkb@conrad.de](mailto:tkb@conrad.de)**  
**Mon to Thurs 8.00am to 4.30pm**  
**Fri                    8.00am to 2.00pm**

## Intended Use

The Body Fat Scales Bettina are used for determining and displaying the body weight and the body fat share of a person. The device was conceived for use in the home.

**The scales are not suitable for medical or therapeutic purposes. These scales do not meet the requirements of the medical products guideline 93/42/EWG.**

Any other use than that described above can lead to damage to the product and, in addition, involves certain risks. No part of the product may be modified or converted in any way. Do not open the housing!

The safety instructions must be observed without fail!

# Table of Contents

Introduction .....	15
Intended use .....	15
Table of contents .....	16
Safety instructions and danger warnings.....	16
Insertion or replacement of battery .....	19
Setting up scales .....	20
Setting the personal data.....	20
Weight and body-fat measurement .....	21
Use as personal scales .....	24
Handling .....	24
Servicing and maintenance .....	24
Disposal .....	25
Elimination of disturbances .....	25
Technical data .....	26
Scope of delivery .....	26

## Safety instruction and danger warnings

The exclamation mark in a triangle indicates that the information provided next to it is of importance in the operating instructions. Read through the operating instructions completely before setting up the device; they include important information for correct operation.

The guarantee will lapse in case of damage caused by failure to comply with these instructions! We shall not be liable for any consequent loss!

We do not accept liability for personal injury or damage to property caused by incorrect handling or non-observance of the safety instructions! Any claim to warranty will lapse in such cases.

For safety and licensing reasons (CE), unauthorised conversion and/or modification of the system is not permitted.

Make sure that the device is set up correctly. Observe these operating instructions in doing so.

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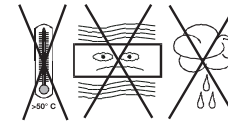
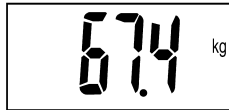
The scales may not be subjected to direct, intensive sunlight, extreme temperatures, strong vibrations, high degrees of humidity, dampness or heavy mechanical stress. Protect the scales from adverse ambient conditions such as, for example, dust or excessive ambient temperatures ( $> +50\text{ }^{\circ}\text{C}$ ) or low ( $< -10\text{ }^{\circ}\text{C}$ ).

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Only a 9-V block battery may be used as a source of power. Never attempt to operate the device with a different source of voltage or with other types of battery.

Do not place the device directly on or next to strong electric or magnetic fields as these could influence the measuring result. Such fields occur in the vicinity of cordless telephones, for example.

The scales are not suited for commercial or professional application, e.g. in fitness studios, hospitals, doctor's offices, etc.

The scales are not suitable for medical or therapeutic purposes.

If you are a diabetic, only use the scales after consulting your physician.

Do not use this product if you have a pacemaker, medicine pump or any other medical implant in or on yourself.

Do not place the scales on slippery, damp or wet floors. There is a risk of the scales sliding away or your slipping when you step on to them.

Do not step too close to the edge of the scales, do not jump on to or off the scales as they could tilt and cause injury.

Do not step on the scales with damp or wet feet. Do not step on the scales when you are wet. There is the risk of slipping and injuring yourself. Dry yourself prior to using the scales.

The manufacturer or supplier do not assume any responsibility for incorrect displays, measurements or the consequences resulting from incorrect displays or measurement.

Handicapped persons and children should not use the scales without proper assistance. There is a risk of injury.

Do not overload the scales. Persons weighing more than 150 kg or 330lb may not use the scales.

The scales may only be used to determine the weight and body fat of persons. The scales are not suited for other weight or body fat measurements, e.g. of animals.

Do not step on to the display as it could be damaged.

The scales are approved for operation in closed rooms only.

The scales are not a toy and therefore have no place in the hands of children.

Do not place any vessels (e.g. flower pots) on the device. Never pour out liquids above the device. In case any liquid penetrates the inside of the scales nevertheless, remove the batteries immediately and contact an expert.

Do not leave packaging material lying around carelessly. In the hands of children, plastic foil and/or bags and polystyrene parts etc. can turn into dangerous toys.

**If you are not certain of how to use the scales correctly or in case of any questions not covered by the operating instructions, please contact our technical support or another expert. Please consult a specialist in case of doubts about correct operation or the product safety.**

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## Inserting the battery / replacing the battery

**Make sure that the polarity is correct when inserting the batteries (pay attention to plus and minus!).**

**Batteries have no place in the hands of children. Do not leave the batteries lying around in the open; there is the risk of their being swallowed by children or domestic animals. Danger to life! If swallowed, contact a doctor immediately.**

**Batteries may not be shorted, opened or thrown into fire.**

**Batteries may not be recharged. There is a risk of explosion!**

**Only use new batteries of the correct size and the recommended type.**

**Chemical liquids can leak from old or used batteries and damage the device. Therefore, remove batteries from the device if you do not require it for a longer period (> 1 month)**

**Leaking or damaged batteries can lead to caustic burns of the skin. Therefore, use suitable protective gloves in this case.**

- The battery compartment is located on the bottom of the scales. Turn the scales over to open it. So as not to scratch the top or the display, the scales should be placed on a soft surface.
- Open the battery compartment. Fasten the battery clip on the battery (9-V block battery) with the correct polarity.
- Now fasten the 9-V block battery to the clamps of the battery compartment cover (see illustration).
- Close the battery compartment again.

It is necessary to replace the batteries when "LO" appears on the display of the scales or if you can no longer turn on the scales.

When you replace the batteries, replace them only with new ones and those of the same type.

## Setting up the scales

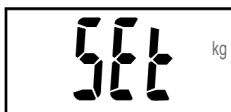
Place the body fat scales horizontally on to a level, hard surface (e.g. tiled floor). Placing the scales on a soft floor, e.g. carpet, or in an incorrect position may falsify the measuring result. Set up the scales in such a way that they stand safely.

Do not place the scales on slippery, damp or wet floors. The scales may slide away when you step on or off them.

## Setting your personal data

In order to perform a body fat analysis, you first of all have to store your personal data. The scales allow you to store the data of up to 12 person.

Press the key . The indication "SET" appears on the display. The indication "kg" or "lb" flashes. Select kilogram (kg) or pound (lb) using the keys  $\triangle$  or  $\nabla$  the desired unit of weight. The selection of the unit of weight, e.g. kg, always applies for all storage places.

Press the key . After that, the current storage place flashes on the display. Select the desired memory slot with the up/down buttons  $\triangle$   $\nabla$ .

In order to enter or alter your personal data, you have to press the button  once again.


The gender display starts to flash. Select your gender with the up/down buttons  $\triangle$   $\nabla$ .

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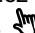
The display for the height of the body (cm or ft in) starts to flash. You can now set the correct body height with the up/down buttons  . Press one of the buttons and keep it depressed in order to start the fast skip mode. Then press the button  to store the value.

The age display starts flash. You can now set the correct age with the up/down buttons  . Press one of the buttons and keep it depressed in order to start the fast skip mode. Then press the button  to store the value.

The keys stop flashing. The scales turn off automatically after approx. 3 seconds.

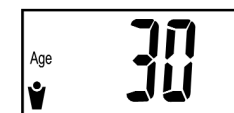If, when entering your personal data, no button is pressed for approx. 30 seconds, the scales will switch off automatically.

If the unit of weight kilogram (kg) is selected, the body height is automatically set to cm.

## Weight and body fat measurement

**Do not use this product if you have a pacemaker or another medical implant in or on yourself.**

**The scales are not suitable for medical or therapeutic purposes. Do not set up any diet or training programmes on the basis of the measuring results. Always consult a physician or another qualified person for this purpose.**


**The measuring result for body fat is subject to many influences. You should observe the following instructions in order to obtain an accurate measuring result.**

Each item of clothing you wear falsifies the measuring result with respect to your weight and body fat.

Do not wear any socks when measuring. Step on the scales with bare feet. The soles of your feet should be clean but not wet.

The soles of the feet have to stand correctly on the electrodes (round metal areas on the weighing platform). The two feet and the thighs may not touch each other.

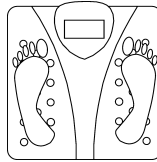The measuring result is influenced strongly after intensive physical activity, during illness, through heavy loss of water (e.g. after heavy consumption of alcohol) after eating and drinking (wait for approx. 2 hours here.), after a bath, in the case of women during menstruation or pregnancy, shortly after rising.

If possible, always take the measurements at the same time of day and under the same conditions.


Press the key. The scales switch  $\triangle$  on  $\nabla$ . The storage place number begins to flash. When the storage place numbers flash, using the up/down buttons  $\triangle$   $\nabla$ , select your personal storage place number.

The display switches to "0.0 kg" (3). After that, the scales are ready for the weight to be measured.

Step on to the scales and stand still until measuring is complete. Weight measurement is complete as soon as the display does not alter anymore.

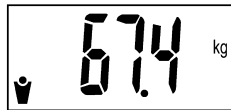The scales automatically change over to measuring the body fat now. Measuring is complete as soon as the display "8888" goes out.

The results flash in the display for approx. 10 seconds in the display. The scales switch off automatically. You can now step off the scales.

You can now compare the values with the body fat table.

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Age	Woman				Man			
	Thin	Normal	Fat	Very fat	Thin	Normal	Fat	Very fat
10-16	<18%	18– 28%	29– 35%	> 35%	<10%	10– 18%	19– 23%	> 23%
17-39	<20%	20– 32%	33– 38%	> 38%	<12%	12– 20%	21– 25%	> 25%
40-55	<23%	23– 35%	36– 41%	> 41%	<13%	13– 21%	22– 26%	> 26%
56- 85	<24%	24– 36%	37– 42%	> 42%	<14%	14– 22%	23– 27%	> 27%

Please read the following information through thoroughly in order to understand the application/informative value of the table correctly.

As the fat and water values determined in percent by the scales are subject to many factors and influences, the displayed values have informative value only and may not be taken as absolute values. You cannot necessarily conclude whether the fat share lies within an ideal range from the physical constitution of a person. Nor do persons of the same age, gender and physique have to have the same fat values. Please consult your physician in order to determine and attain the best fat value for yourself. On the basis of the values determined, the scales can help you to monitor the tendency of change over a longer period of time.


The values in the table have no informative value for the following persons:

- Children below the age of 10 or persons above the age of 85
- Persons suffering from fever, oedema and osteoporosis
- Persons who take medication for the heart
- Dialysis patients
- Pregnant women
- Professional athletes who train for more than 10 hours a week

## Use as personal scales

You can also use the body fat analysis scale as a normal set of scales.


Press the  key briefly. The scales switch on. The indication "8888" appears for approx. 2 seconds (do not step on the scales yet).

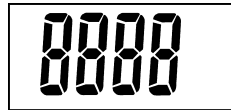As soon as "0.0 kg" or "0.0 lb" appears in the display, the scales are ready for the weight measurement. Step on to the scales and stand still until weight measuring is complete.

Weight measurement is complete as soon as the display does not alter any longer and then flashes 3 x. You can now step off the scales. The weight display remains for approx. The scales then turn off automatically. (Illustration 6)
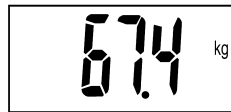
## Handling

Never switch on the device immediately after it has been taken from a cold to a warm environment. The condensation water which develops may destroy the product.

**Allow the device to reach room temperature before switching it on. Wait until the condensation has evaporated.**

Avoid using the device under unfavourable ambient conditions. The following are unfavourable ambient conditions:

- excessive humidity (> 90 % rel. humidity), moisture, dust, combustible gases, vapours or solvents
- excessive ambient temperatures (> +40 °C), ambient temperature too low (< +5 °C)

## Servicing and maintenance

Observe the following safety instructions before the device is cleaned or maintained:

Only a specialist who is familiar with the risks and relevant regulations may carry out repairs or open the device.

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The device is maintenance-free. For this reason, you should never open it. If you do open it, the warranty will lapse.

Only clean the outside of the body fat analysis scales with a soft, slightly damp cloth or a dry brush. Do not use any aggressive cleaning agents or chemical solutions as these may damage the device or impair its function.

## Disposal

Dispose of the used body fat analysis scales in accordance with the valid statutory regulations. Remove the batteries from the scales beforehand.

example

### Disposal of flat batteries / accumulators

You, as the ultimate customer, are legally obliged (**Regulation regarding Flat Batteries**) to return all used batteries and accumulators; **disposal in the household waste is prohibited!**

Batteries / accumulators which contain hazardous substances are marked with the symbols on the side. These symbols indicate that it is prohibited to dispose of these batteries in the household waste. The abbreviations for the decisive heavy metal are: **Cd** = cadmium, **Hg** = mercury, **Pb** = lead. You can hand in your used batteries/accumulators at the official collection points of your community at no cost or everywhere where batteries/accumulators are sold.

**You thus fulfil the legal requirements and make your contribution to the protection of the environment!**

## Troubleshooting

By purchasing this system you have acquired a product which has been designed to the state of the art and is operationally reliable.

However, problems and malfunctions may occur nevertheless.

For this reason, we wish to describe to you here how you can eliminate possible disturbances yourself.

**Always observe the safety instructions!**

Problem	Cause	Solution
Indication "FFFF" in the display	Measuring range has been exceeded (over 150kg or 330lb)	No measurement is possible in this case
Indication "LO" appears in the display	The battery voltage is too low.	Replace the battery
Indication "En" appears in the display	Socks or shoes were worn during measurement	Only measure the body fat with bare feet
	The weight and fat values determined can not match the values entered on the storage place.	Check the values entered for your personal storage place. Was the correct storage place selected?

Repairs other than those just described should only be performed by an authorised specialist.

## Technical data

Batteries:	9-V block battery (6AM6)
Power consumption:	6mA
Measuring range weight / division:	2kg (min) to 150kg (max.) / 0.1kg 4.4lb (min) to 330lb (max.) / 0.2lb
Measuring range - body fat / division:	3% - 75% / 0,1%
Adjustable range - body height:	100cm - 250cm
Adjustable age:	5 - 85

## Scope of delivery

Body fat scales  
9-B block battery  
Operating instructions

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© Copyright 2004 par Conrad Electronic GmbH. Imprimé en Allemagne.